

**SARS-Related Coronavirus 2 (SARS-CoV-2)  
Lineage JG.3; Omicron Variant  
Isolate: USA/NY-Wadsworth-23067147-01/2023  
Culture Fluid (Heat Inactivated) (1 mL)**

## PRODUCT DESCRIPTION:

SARS-Related Coronavirus 2 Lineage JG.3; Omicron Variant (Isolate: USA/NY-Wadsworth-23067147-01/2023) is an enveloped, positive-sense single-stranded RNA virus from the *Coronaviridae* family and the *Betacoronaviridae* genus.

SARS-CoV-2 isolate USA/NY-Wadsworth-23067147-01/2023 was collected in November 2023 in New York, USA. Under the nomenclature system introduced by GISAID (Global Initiative on Sharing All Influenza Data), SARS-CoV-2, isolate USA/NY-Wadsworth-23067147-01/2023 is assigned lineage JG.3. The complete genome of USA/NY-Wadsworth-23067147-01/2023 has been sequenced (GISAID: EPI\_ISL\_18699678).

Each frozen aliquot contains 1 mL of heat inactivated viral culture fluid. The pre-inactivation titer was determined from an infectious aliquot.

Viral inactivation is verified after heat inactivation by the absence of viral growth in tissue culture-based infectivity assays.

## INTENDED USE:

Heat inactivated viral culture fluids are sold as consumable testing materials. The suitability and performance characteristics should be determined by your laboratory for each intended usage.

These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under section 351 of the Public Health Service Act or for any other product intended for administration to humans.

## FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

## BIOSAFETY:

Please consult your institution's regulations regarding the use of this product. For a detailed discussion on biological safety see the current edition of Biosafety in Microbiological and Biomedical Laboratories (BMBL), published by the CDC.





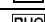


## PRECAUTIONS:

- Use Universal Precautions, this product is **potentially biohazardous**.
- Repetitive freezing and thawing is not recommended (aliquot material if necessary).
- To avoid cross-contamination, use separate pipette tips for all reagents.

## RECOMMENDED STORAGE:

Heat inactivated viral culture fluids should be stored at -65°C or below.

PI0810696CFHI  
Revision: 00  
Effective Date: 04/29/2024

 REF	Catalog Number		Temperature Limitation
 LOT	Batch Code		Expiration Date
 RUO	For Research Use Only		Biological Risk
	Manufacturer		

www.ZeptoMetrix.com

ZeptoMetrix LLC • 878 Main Street, Buffalo, NY 14202 USA • Tel (800) 274-5487

This product was manufactured in a facility which has a Quality Management System that is ISO 13485 certified.

PCA# N/A  
Page 1 of 1