

**SARS-Related Coronavirus 2 (SARS-CoV-2)
Lineage B.1.1.529; Omicron Variant¹
Isolate: USA/CO-CDPHE-2102544747/2021
Culture Fluid (Heat Inactivated) (1 mL)**

Catalog Number:
0810643BCFHI

PRODUCT DESCRIPTION:

Severe acute respiratory syndrome-related coronavirus 2 (SARS-CoV-2), isolate hCoV-19/USA/CO-CDPHE-2102544747/2021 was isolated from a human on December 29, 2021, in Colorado, USA. Under the nomenclature system introduced by GISAID (Global Initiative on Sharing All Influenza Data), SARS-CoV-2, isolate hCoV-19/USA/CO-CDPHE-2102544747/2021 is assigned lineage BA.2 (Pango v.3.1.20 2022-02-28), Omicron (BA.2-like) (Scorpio) and GISAID clade GRA using Phylogenetic Assignment of Named Global Outbreak lineages (PANGO) tool. The complete genome of the clinical isolate of SARS-CoV-2, hCoV-19/USA/CO-CDPHE-2102544747/2021 has been sequenced (GISAID: EPI_ISL_8643930).

Each frozen aliquot contains 1 mL of heat inactivated viral culture fluid. The pre-inactivation titer was determined from an infectious aliquot.

Viral inactivation is verified after heat inactivation by the absence of viral growth in tissue culture-based infectivity assays.

INTENDED USE:

Heat inactivated viral culture fluids are sold as consumable testing materials. The suitability and performance characteristics should be determined by your laboratory for each intended usage.

These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under section 351 of the Public Health Service Act or for any other product intended for administration to humans.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

BIOSAFETY:

Please consult your institution's regulations regarding the use of this product. For a detailed discussion on biological safety see the current edition of Biosafety in Microbiological and Biomedical Laboratories (BMBL), published by the CDC.

PRECAUTIONS:





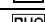


- Use Universal Precautions, this product is **potentially biohazardous**.
- Repetitive freezing and thawing is not recommended (aliquot material if necessary).
- To avoid cross-contamination, use separate pipette tips for all reagents.

RECOMMENDED STORAGE:

Heat inactivated viral culture fluids should be stored at -65°C or below.

¹The following reagent was deposited by the Centers for Disease Control and Prevention and obtained through BEI Resources, NIAID, NIH: SARS-Related Coronavirus 2, Isolate hCoV-19/USA/CO-CDPHE-2102544747/2021 (Lineage B.1.1.529, BA.2; Omicron Variant), NR-56520.

PI0810643BCFHI
Revision: 01
Effective Date: 10/17/2023

	Catalog Number		Temperature Limitation
	Batch Code		Expiration Date
	For Research Use Only		Biological Risk
	Manufacturer		

www.ZeptoMetrix.com

ZeptoMetrix LLC • 878 Main Street, Buffalo, NY 14202 USA • Tel (800) 274-5487

This product was manufactured in a facility which has a Quality Management System that is ISO 13485 certified.

PCA# 23-351
Page 1 of 1